


24H ORAL HYGIENE KIT

B-Care Oral Care

An Innovative Approach to Oral Care

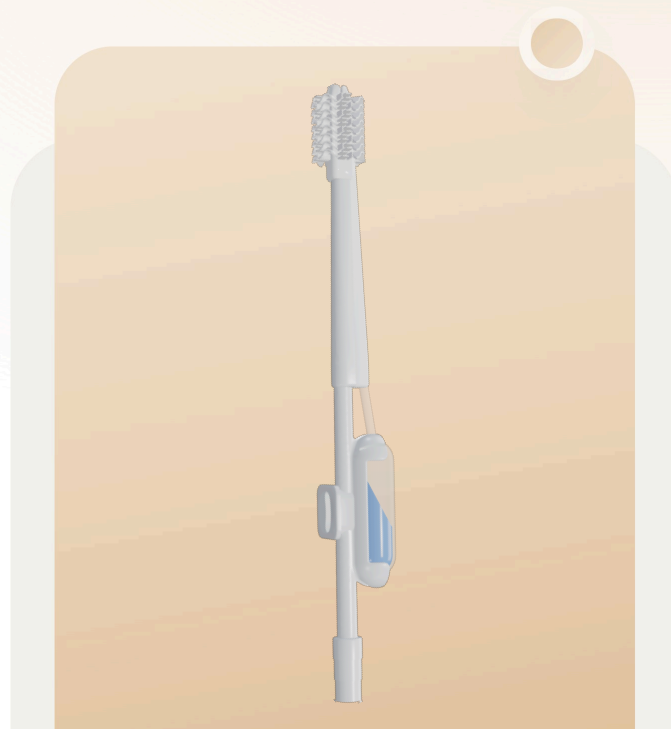

• MECHANISM

A soft 360 degree brush, a 4 ml pipette of mouthwash in the handle and an integrated suction channel: the mouthwash goes straight to the brush head, and the suction takes it away, all in one hand.

Brush	360 degree, soft TPE
Brush head	16 × 18 mm
Fluid	4 ml pipette, mouthwash
Dosing	Thumb on the pipette
Suction	Integrated, thumb control
Connection	Medical suction unit
Length	178 mm
Patients	Intubated adults

• APPLICATION

- 1 Fill**
Fill the 4 ml pipette with mouthwash.
- 2 Insert**
Place it in the cradle, the tip in the fluid inlet.
- 3 Connect**
Attach the handle to the medical suction unit.
- 4 Clean**
Press the pipette with the thumb and gently brush.

Every 8 h
one cleaning

3 packs
one blister, 24 hours

• BENEFITS

-  **360° soft, rubber-like (TPE) brush:** easy access to the patient's mouth
-  **Disinfectant-moisturizing fluid** stored in handle, dispensed via brush tip
-  **Simple tracking of 24-hour oral hygiene protocol**
-  **Single-hand fluid usage**

Toothbrushing is an essential practice for ventilated adults (SHEA/IDSA/APIC 2022).

• ORDERING DATA

PRODUCT
B-Care Oral Care
24h Oral Hygiene Kit

PRODUCT NUMBER
BC241

UNIT
1 blister (set of 3)
Size/unit 230 × 200 mm

OUTER CARTON
205 × 320 × 235 mm
20 blisters/carton

SET
B-Care Oral Care
Toothbrush
two pipettes (4 ml),
swab, saliva ejector

MANUFACTURER
Airway Medix S.A.
Szara 10, 00-420
Warsaw, Poland

STORAGE
+5 to +30 °C
shelf life 2 years

CONFORMITY
CE, EU MDR 2017/745

ORDERS
sales@airway-medix.com
airway-medix.com



Product page
IFU, tutorial, downloads