

24H ORAL HYGIENE KIT


B-Care Oral Care

An Innovative Approach to Oral Care

• ONE HAND

**Single-hand fluid
usage.**

Disinfectant-moisturizing fluid
stored in handle, dispensed via
brush tip

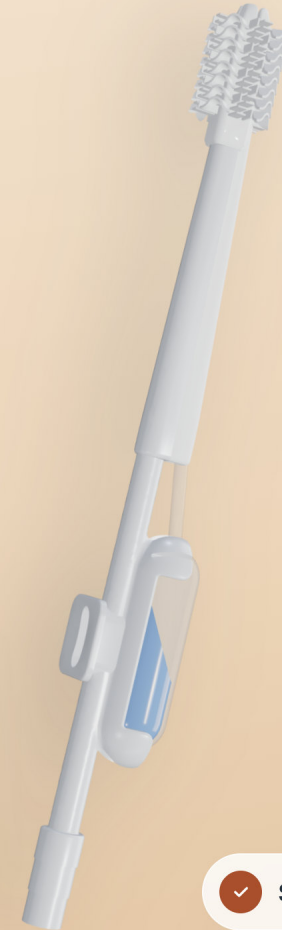


✓ Single use · intubated adults

24H ORAL HYGIENE KIT

B-Care Oral Care

● KEY DATA

Brush	360 degree, soft TPE
Fluid	4 ml pipette
Dosing	Thumb on the pipette
Suction	Integrated
Pack	3 per day, 1 every 8 h
Set	Toothbrush, 2 pipettes (4 ml), swab, saliva ejector
Use	Single use, adults

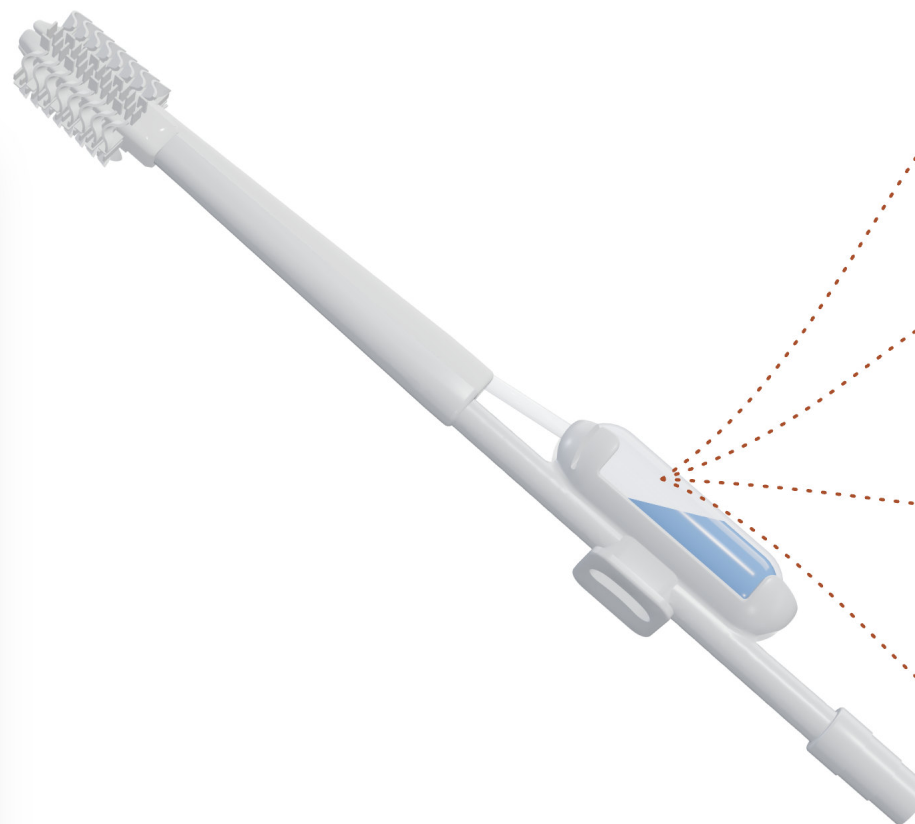

360° soft brush
360° soft, rubber-like (TPE) brush: easy access to the patient's mouth

Fluid in the handle
Disinfectant-moisturizing fluid stored in handle, dispensed via brush tip

24-hour protocol
Simple tracking of 24-hour oral hygiene protocol

One hand
Single-hand fluid usage

TECHNICAL SPECIFICATION

Small. Simple. Single use.


16 × 18 mm brush head



4 ml pipette



Suction connection

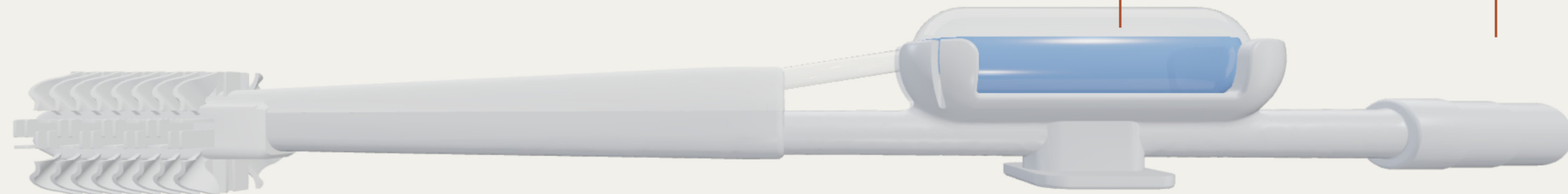

178 mm overall length

● DEVICE

Brush	Cylindrical, 360 degree
Bristles	Soft TPE, atraumatic
Brush head	16 × 18 mm
Length	178 mm
Handle	PP, suction channel
Build	2 parts, bonded

● FLUID AND SUCTION


Pipette	4 ml, replaceable
Size	12.4 × 40 mm, PE
Fluid	Mouthwash
Dosing	Thumb on the pipette
Suction	Integrated, thumb control
Connection	Medical suction unit

● SUPPLY AND REGULATORY

Supply	1 blister (set of 3), 24 hours
Interval	Every 8 hours
Patients	Intubated adults
Use	Single use, non-sterile
Shelf life	2 years, +5 to +30 °C
Conformity	CE, EU MDR 2017/745

B-CARE ORAL CARE
AIRWAY MEDIX S.A.

EXPLODED VIEW
2 PARTS + PIPETTE · NOT TO SCALE

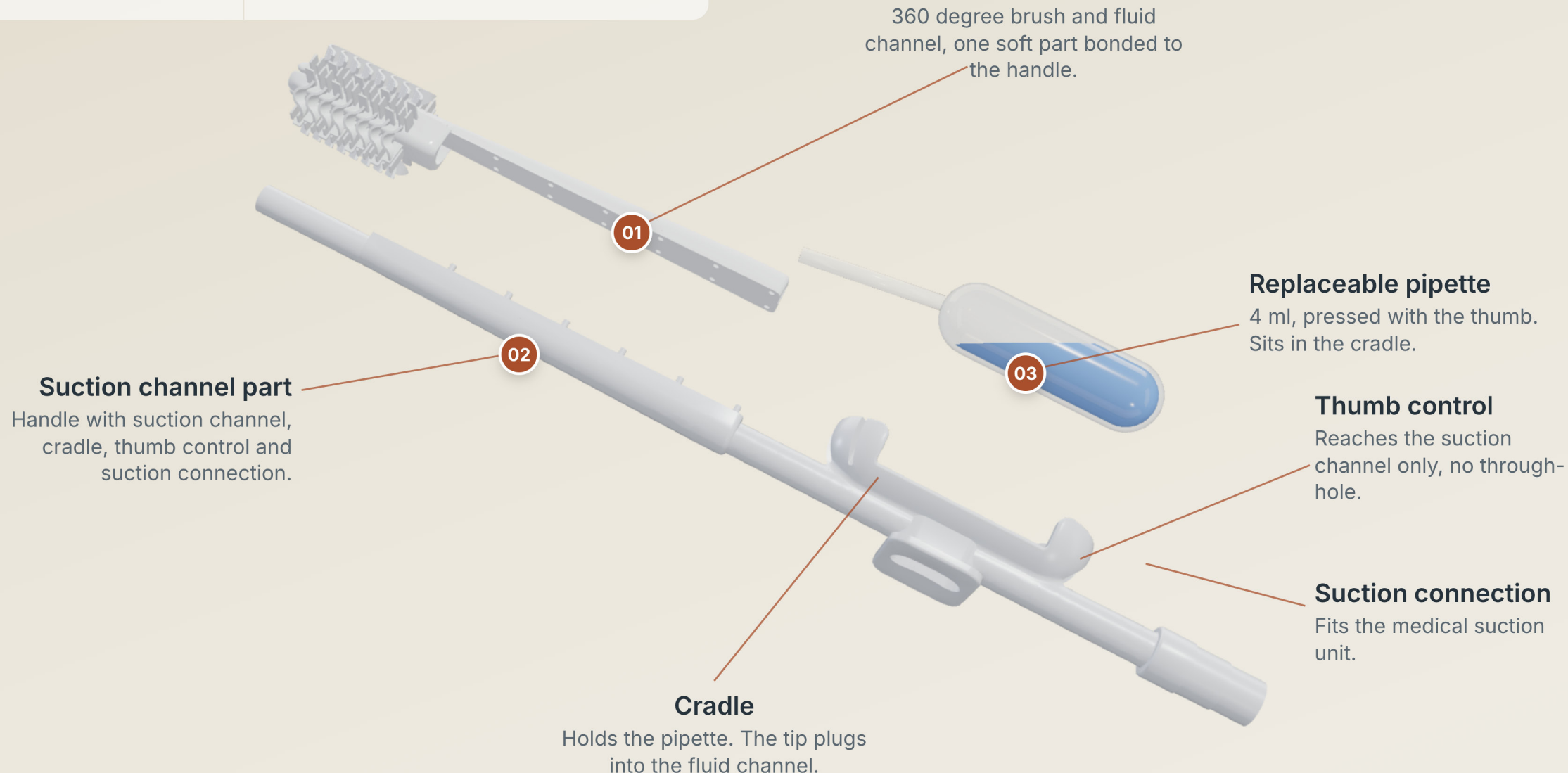

IN USE

Ready in four steps.

**Every 8 h**

one cleaning

3 packs

one blister, 24 hours

1**Fill**

Fill the 4 ml pipette with mouthwash.

2**Insert**

Place it in the cradle, the tip in the fluid inlet.

3**Connect**

Attach the handle to the medical suction unit.

4**Clean**

Press the pipette with the thumb and gently brush.

Toothbrushing is an essential practice for ventilated adults (SHEA/IDSA/APIC 2022).

COMPARED

Only B-Care has all four.

FEATURE	B-CARE	FIVE LEADING KITS
Fluid on board Pipette sits in the handle.	✓	None
360 degree brush Cleans all around.	✓	None
Flexible suction In the same hand.	✓	None
Compact kit One pack per shift.	✓	2 of 5



One dedicated partner per market.

We work with one distributor in each market. You talk to us, you sell, we support you.



1

Talk to us.

Tell us about your market and the hospitals you serve.

2

Agree on the terms.

Volumes and conditions for your market, set together. One partner per market: in yours, that is you.

3

Get your support material.

Registration documents, IFU, tutorials, trailer and one-pagers, ready to use.

Let's simplify oral care in your country.

We work with one distributor in each market. You talk to us, you sell, we support you.

sales@airway-medix.com

airway-medix.com

REVOLUTIONIZING MEDICINE, SIMPLIFYING CARE